

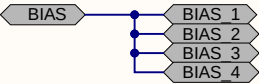
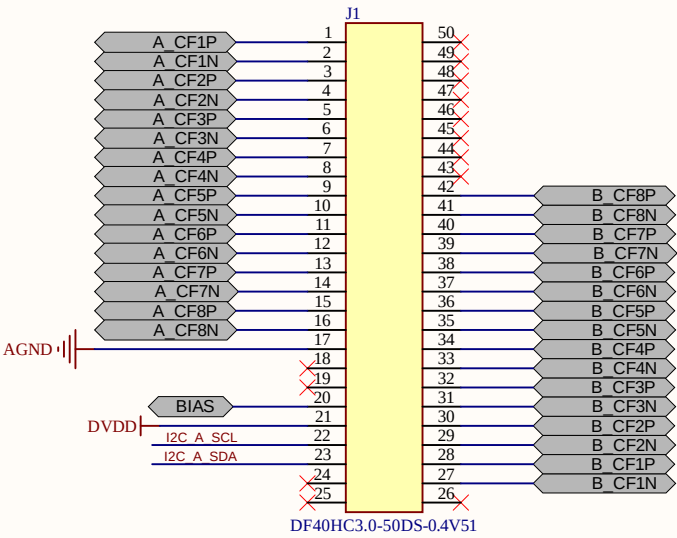


Fully differential configuration

50-pin board-to-board connector (male)

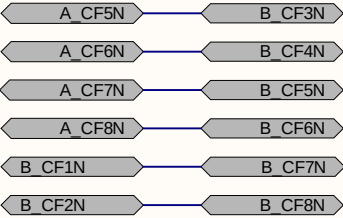

40-pin board-to-board connector (male)

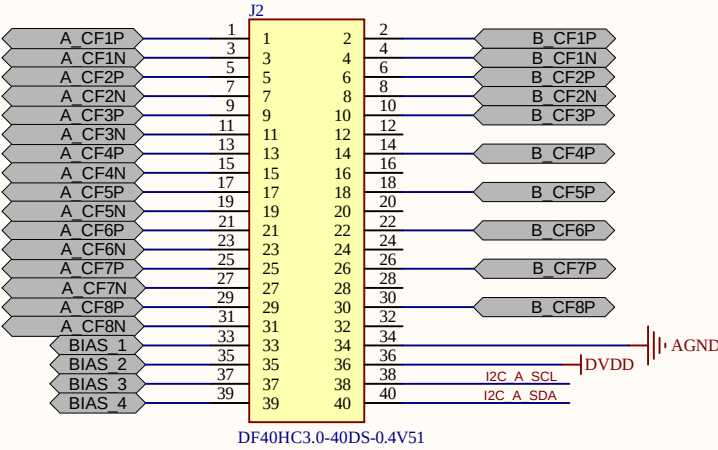

- For single ended configuration we need:
- Common reference (16 channels max)
 - AGND which is used to set the body bias in case of bipolar supply
 - BIAS in case of unipolar supply and pullups
 - AVDD and AVSS (in case of active electrodes)
 - Two shared GPIOs for future use



Eidgenössische Technische Hochschule Zürich
Swiss Federal Institute of Technology Zurich

Project:

BioGAP EMG shield

Drawing number: *

Rev: v2.0

Format:

Laboratory: *

Sheet: 01_channel_share_config.SchDoc

Date: 03.07.2025 17:19:25

A3

Drawn by: *

Page * of *

File: C:\Users\ls-fre\Github\emg_sleeve\Hardware\EMG_shield_channel_share\01_channel_share_config.SchDoc